

STUDY PUBLICATION OF SELECTED RESULTS

EFFECTIVENESS AND LONG-TERM EFFECT OF A KNEE SUPPORT IN PATIENTS WITH CHRONIC INSTABILITY FOLLOWING ACL RECONSTRUCTION

BACKGROUND

The surgical reconstruction of a torn cruciate ligament with subsequent rehabilitation is the most commonly pursued treatment in young, active patients. [1] The annual incidence rate for ACL reconstructions per 100,000 is 68.6 in the US [2], 58.2 in New Zealand [3], 52.0 in Australia [4] and 32.0 in Sweden. [5] The incidence is higher in populations playing team sports than for people who do not exercise. [6]

In the post-operative setting, reports show medium-term to long-term impairment of and restriction in knee functionality, in addition to the risk of rerupture [7] or the development of osteoarthritis of the knee [8].

Clinical studies suggest that potential consequences following cruciate ligament reconstruction may include persistent thigh muscle deficits [9,10], abnormal gait [11] and reduced physical activity [12]. Activity levels also seem to decrease over time [13], with excessive body weight playing a part as well [14].

A long-term decrease in the quality of life related to knee function, fear of recurring cruciate ligament tears and diminishing confidence are potential developments that have been observed in patients following cruciate ligament ruptures [10,15–17]. Studies report on the use of elastic supports in individual patients during the rehabilitation phase and when returning to exercise following ACL reconstruction [18]. It has been discussed that supports may improve/normalize gait by improving proprioception and sensorimotor control [19,20], thus increasing knee joint function and the patient's confidence in their own knee [21,22].

The objective of the study was to examine the stabilizing effect of the GenuTrain knee support in patients with chronic instability (at least 5 months after surgery) following ACL rupture and ligament reconstruction at baseline as well as 6 weeks after wearing the product.

SOURCE

Sole, G., Lamb, P., Pataky, T., Pathak, A., Klima, S., Navarre, P., and Hammer, N. Immediate and six-week effects of wearing a knee sleeve following anterior cruciate ligament reconstruction on knee kinematics and kinetics: a cross-over laboratory and randomised clinical trial. *BMC Musculoskeletal Disorders* (2022) 23:560

STUDY DESIGN

Crossover design for the acute effect; two-arm, randomized, controlled clinical study with 6-week long-term follow-up (evidence level 1b)

METHODOLOGY

Sample: **n=34 patients;**
Part 1: acute effect: n=34 (crossover, randomized)
Part 2; 6 weeks of wearing the product:
n=17 with support = BG = intervention group,
n=17 without support = KG = control group

Intervention group:
Age: 27 ±7 years,
Height: 173.0 ±10 cm,
Weight: 72.9 ±10.7 kg,
BMI: 24.4 ±3.2, sex: m:w=10:7,
Time since surgery: 15 months (7–44)
Wearing duration of the support/day:
at least one hour
*Tegner Activity Scale before the injury: 7 (3–10)
*Tegner Activity Scale since the injury,
after surgery: 4 (2–9)

Control group:
Age: 26 ±7 years,
Height: 173.0 ±01 cm,
Weight: 80.4 ±11.1 kg,
BMI: 26.7 ±2.4, sex: m:w=8:9;
Time since surgery: 16 months (6–53)
*Tegner Activity Scale before the injury: 9 (6–10)
*Tegner Activity Scale since the injury,
after surgery: 5 (3–9)

*Tegner Score: 0 – Unable to participate in sports/activities because of knee problems; 5 – Recreational exercise is possible, jogging on uneven surfaces 2x week, 10 – Competitive national or international sports, team sports (football, rugby)

Sole, G., Lamb, P., Pataky, T., Pathak, A., Klima, S., Navarre, P., and Hammer, N. Immediate and six-week effects of wearing a knee sleeve following anterior cruciate ligament reconstruction: a cross-over laboratory and randomised clinical trial. *BMC Musculoskeletal Disorders* (2021) 22:655

METHODOLOGY

Test supports:	GenuTrain® (Bauerfeind AG)
Measurement systems and test method:	<p>Horizontal jump: Single-leg jump from standing on one leg and landing on one leg, average of 3 attempts, healthy side, injured side without and with the support</p> <p>Step-down hop test: Single-leg jump from a box 30 cm in height onto a load sensing platform positioned 15 cm in front. Landing on one leg, immediately followed by a forward jump from the load sensing platform to a distance of 60 to 70% of the maximum horizontal jumping length. Average of 3 attempts, healthy side, injured side without and with the support.</p> <p>LSI Lateral Symmetry Index (injured side/healthy side * 100)</p> <p>IKDC-SKF International Knee Documentation Classification</p> <p>Subjective Knee Form = function, max. points: 100</p>
Study period:	<p>1st measurement: Acute effect, measurement on the day the support was handed out followed by 6 weeks of wearing the support for BG, n=17 and KG n=17 without a support</p> <p>2nd measurement: 6 weeks after 1st measurement</p>
Inclusion criteria:	<p>Patients with ACL rupture and ACL reconstruction at least 5 months to 5 years in the past, functional deficits, measured using the IKDC-SKF (Intern. Knee Documentation Committee Subjective Knee Form); 40–80 of 100 points</p>
Exclusion criteria:	<p>Patients with revision procedures or previous ACL ruptures on the other knee</p> <p>Patients with problems related to the pelvis or lower back as well as the lower extremities</p> <p>Conditions that required medical treatment during the past 6 months or that caused a restriction in daily activity</p> <p>Neurological or cardiovascular conditions</p> <p>BMI above 30 or IKDC-SKF value < 40 or > 80</p>

RESULTS

Single-leg horizontal jump, acute effect:

During jumps with the support on the injured leg, the distance increased significantly by 3.6% (95% CI 0.4–6.8%, $p=0.025$) compared with jumping without a support on the injured leg.

A reduction in different jumping abilities was also observed between the healthy and the injured side of -9.3% (-12.4%, -6.1%) without a support to -6.0% (-9.2%, -2.8%) when wearing the support. During the acute phase, the deficit on the injured side compared with the healthy side decreases by a third when wearing the GenuTrain. This corresponds to an increase in jumping ability of 5 cm for the injured leg when a support is being worn. (Fig. 1)

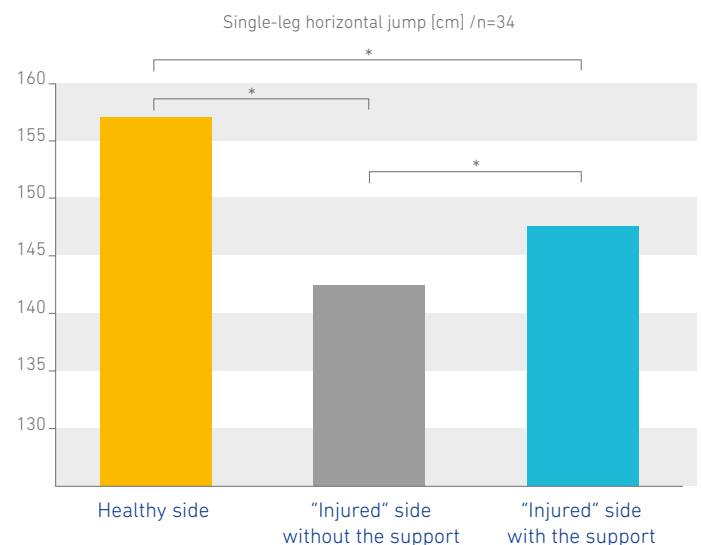


Fig. 1: Single-leg horizontal jump Y axis = jumping length [cm], ($\alpha<0.05$; power, $\beta=80\%$; one-way ANOVA)

SINGLE-LEG HORIZONTAL JUMP, LONG-TERM EFFECT:

After 6 weeks of intervention, no significant further increase in jumping ability was measured during the single-leg jump using the injured side when comparing "without" and "with" support. However, a trend can be noticed where the symmetry index in the intervention group (jumping ability; jumping distance with the injured side compared with the healthy side) increased from 92.0 per cent to 95.2 per cent. (Fig. 2)

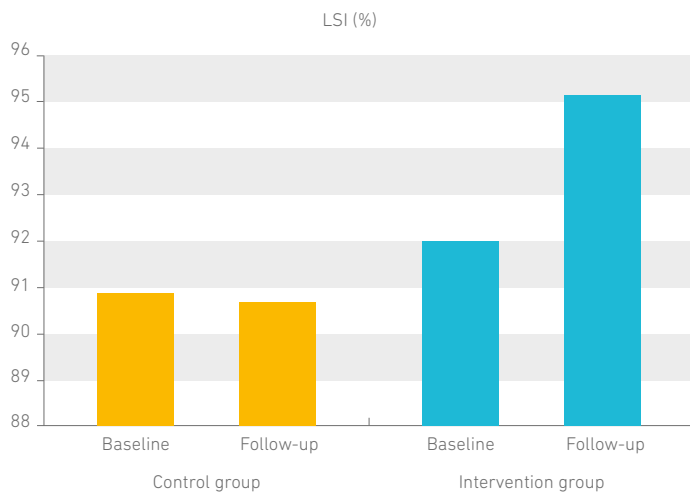


Fig. 2: LSI (Lateral Symmetry Index), comparing the healthy side with the injured side during the acute phase (baseline) and after 6 weeks of wearing the support (follow-up); n=12 KG; n=12 BG

IKDC SCORE, KNEE FUNCTION SELF-ASSESSMENT QUESTIONNAIRE:

There were no significant differences in the overall IKDC assessment between the control group and the intervention group. During the detailed review, however, a significant difference was observed between the control group and the intervention group when asked about an impaired function based on a "locked knee" and/or the knee popping and snapping. In the control group, the number of events was unchanged, whereas the intervention group showed a decrease in the event rate from 71% to 27% after 6 weeks of wearing the support.

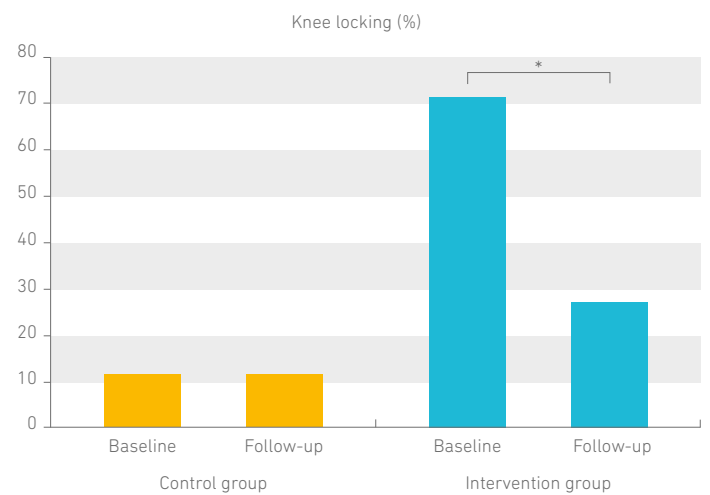


Fig. 3: Comparing knee function before being given the support (baseline) and knee function after 6 weeks of wearing the support (intervention group vs. 6 weeks without a support (control group)), ($\alpha < 0.05$; power, $\beta = 80\%$; McNemar's test), n=16 KG; n=15 BG

SINGLE-LEG HORIZONTAL JUMP (STEP-DOWN HOP TEST); ACUTE EFFECT:

During the single-leg jump without the support on the injured leg, the max. knee angle of the injured side is only 91.9% of the max. knee angle of the healthy side. During jumps with the support on the injured leg, the maximum knee angle increased significantly by 6.8% compared with jumping without the support on the injured leg. When jumping with the support on the injured leg, a max. knee angle was recorded that corresponded to 98.2% of the healthy side.

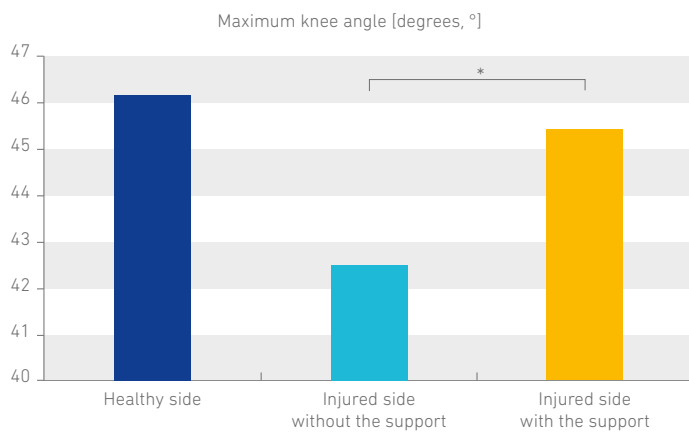


Fig. 4: Step-down hop test, Y axis = max. knee angle [degrees °], ($\alpha < 0.05$; power, $\beta = 80\%$; one way repeated ANOVA)

The max. moments of force during landing were also significantly increased by 6.9% when wearing the support on the injured leg, compared with the injured leg without the support, fig. above.

SINGLE-LEG HORIZONTAL JUMP (STEP-DOWN HOP TEST); LONG-TERM EFFECT:

After six weeks of wearing the GenuTrain, test subjects recorded a reduced standing phase during the step-down hop test with the injured leg (minus 22%). Test subjects without the support recorded no significant change to the injured leg. Results indicate that the GenuTrain can improve jumping performance, over an extended period as well.

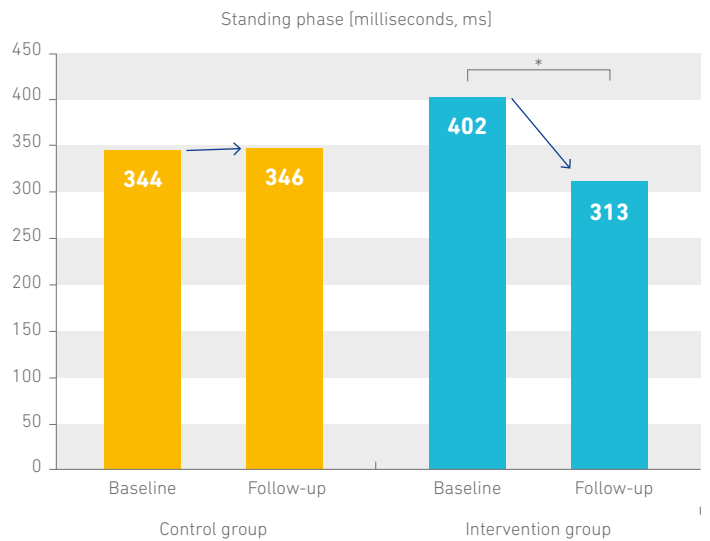


Fig. 5: Step-down hop test, Y axis = standing time [milliseconds, ms], baseline = T0, at the beginning of the test; follow-up = T1 after 6 weeks; ($\alpha < 0.05$; power, $\beta = 80\%$; one way repeated ANOVA)

DISCUSSION

The single-leg horizontal jump is often used to assess recovery in patients following cruciate ligament reconstruction. It is also used in combination with further tests to determine the timing of a safe return to sporting activities [23,24].

Patients in this study demonstrated lateral symmetry (LSI) of LSI $\geq 90\%$ in relation to jumping distance, which is often used as a benchmark for returning to exercise. However, the jumping distance of patients using the injured leg without wearing a support was significantly less in this study than the average distance of 187 cm reported in other studies for patients with comparable injuries and age distribution [25].

On the one hand, wearing the knee support can increase performance by 3.6% during the acute phase. This points to an improvement in the functional status of the "injured" knee.

On the other hand, quite a wide dispersion in the individual data (standard deviation; SEM 3%) is documented. In addition, patients reported wearing the support for very different daily durations over the 6 weeks. The time ranges from 1–8 hours per day.

The wide range of the Tegner Activity Score from 2–9 out of 10 also shows a certain heterogeneity of the groups examined, with some subjects returning to competitive sports (9/10), while others were only just able to walk on uneven surfaces (2/10).

These may be reasons influencing the comparison between the two groups, resulting in further increases in performance not being observed for all test parameters over the total of 6 weeks for all patients in the intervention group. One effect that did continue for the 6 weeks was reflected in the improved symmetry of the injured side compared with the healthy side as well as the notably lower rate of impaired knee function after 6 weeks of wearing the knee support. Furthermore, a significant reduction in the standing phase was measured during the "step-down hop test" which, as a result of the effects, would indicate an improvement of the knee function triggered by the GenuTrain.

Knee supports can be a useful addition to the treatment regimen during rehabilitation and for the ongoing improvement of physical activity.

Based on the evidence of this study data and current existing studies [18–22] examining the effectiveness of knee supports, their use should be determined based on individual assessments and the patient's current recovery progress.

CONCLUSIONS



GenuTrain increases performance and knee coordination



After 6 weeks of wearing the GenuTrain, there was no decrease in the acute effect



After 6 weeks, there were fewer instances of impaired knee function (knee locking) when using GenuTrain



After six weeks, there was a reduced standing phase during the step-down hop test.

REFERENCES

- [1] Rahardja R, Zhu M, Love H, Clatworthy MG, Monk AP, Young SW: Effect of graft choice on revision and contralateral anterior cruciate ligament reconstruction: Results from the New Zealand ACL registry. *The American Journal of Sports Medicine* 2019, 48:63-69.
- [2] Sanders TL, Maradit Kremers H, Bryan AJ, Larson DR, Dahm DL, Levy BA, Stuart MJ, Krych AJ: Incidence of anterior cruciate ligament tears and reconstruction: a 21-year population-based study. *American Journal of Sports Medicine* 2016, 44:1502-1507.
- [3] Sutherland K, Clatworthy M, Fulcher M, Chang K, Young SW: Marked increase in the incidence of anterior cruciate ligament reconstructions in young females in New Zealand. *Australian New Zealand Journal of Surgery* 2019, 89:1151-1155.
- [4] Janssen KW, Orchard JW, Driscoll TR, van Mechelen W: High incidence and costs for anterior cruciate ligament reconstructions performed in Australia from 2003–2004 to 2007–2008: time for an anterior cruciate ligament register by Scandinavian model? *Scandinavian Journal of Medicine and Science in Sports* 2012, 22:495-501.
- [5] Granan LP, Forssblad M, Lind M, Engebretsen L: The Scandinavian ACL registries 2004–2007: baseline epidemiology. *Acta Orthopaedica* 2009, 80:563-567.
- [6] Moses B, Orchard J, Orchard J: Systematic review: Annual incidence of ACL injury and surgery in various populations. *Research in Sports Medicine* 2012, 20:157-179.
- [7] Paterno MV, Rauh MJ, Schmitt LC, Ford KR, Hewett TE: Incidence of second ACL injuries 2 years after primary ACL reconstruction and return to sport. *American Journal of Sports Medicine* 2014, 42:1567-1573.
- [8] Barenius B, Ponzer S, Shalabi A, Bujak R, Norlén L, Eriksson K: Increased risk of osteoarthritis after anterior cruciate ligament reconstruction: a 14-year follow-up study of a randomized controlled trial. *American Journal of Sports Medicine* 2014, 42:1049-1057.
- [9] Petersen W, Taheri P, Forkel P, Zantop T: Return to play following ACL reconstruction: a systematic review about strength deficits. *Archives of Orthopaedic and Trauma Surgery* 2014, 134:1417-1428.
- [10] Kaur M, Ribeiro DC, Lamb P, Webster K, Sole G: Low knee-related quality of life and persistent physical asymmetries in participants up to 10 years post-ACL reconstruction – a cross-sectional study. *Physical Therapy in Sport* 2020, In press: Accepted 11/12/20.
- [11] Kaur M, Ribeiro DC, Theis JC, Webster KE, Sole G: Movement patterns of the knee during gait following ACL reconstruction: a systematic review and meta-analysis. *Sports Medicine* 2016, 46:1869-1895.
- [12] Bell DR, Pfeiffer KA, Cadmus-Bertram LA, Trigsted SM, Kelly A, Post EG, Hart JM, Cook DB, Dunn WR, Kuenze C: Objectively measured physical activity in patients after anterior cruciate ligament reconstruction. *American Journal of Sports Medicine* 2017, 45:1893-1900.
- [13] Spindler KP, Huston LJ, Chagin KM, Kattan MW, Reinke EK, Amendola A, Andrich JT, Brophy RH, Cox CL, Dunn WR et al: Ten-year outcomes and risk factors after anterior cruciate ligament reconstruction: A MOON Longitudinal Prospective Cohort Study. *American Journal of Sports Medicine* 2018, 46:815-825.
- [14] de Oliveira FCL, Roy J-S, Pappas E: ACL injury, physical activity, and overweight/obesity: a vicious cycle? *Knee Surgery, Sports Traumatology, Arthroscopy* 2019, 28:667-669.
- [15] Filbay SR, Ackerman IN, Dhupelia S, Arden NK, Crossley KM: Quality of life in symptomatic individuals after anterior cruciate ligament reconstruction, with and without radiographic knee osteoarthritis. *Journal of Orthopaedic and Sports Physical Therapy* 2018, 48:398-408.
- [16] Filbay SR, Ackerman IN, Russell TG, Macri EM, Crossley KM: Health-Related quality of life after anterior cruciate ligament reconstruction: a systematic review. *American Journal of Sports Medicine* 2013, 42:1247-1255.
- [17] Ardern CL, Taylor NF, Feller JA, Webster KE: Fear of re-injury in people who have returned to sport following anterior cruciate ligament reconstruction surgery. *Journal of Science and Medicine in Sport* 2012, 15:488-495.
- [18] Kaur M, Ribeiro DC, Theis JC, Webster KE, Sole G: Individuals' experiences of the consequences of anterior cruciate ligament reconstruction surgery. *New Zealand Journal of Physiotherapy* 2019, 48:76-93.
- [19] Van Tiggelen D, Coorevits P, Witvrouw E: The use of a neoprene knee sleeve to compensate the deficit in knee joint position sense caused by muscle fatigue. *Scandinavian Journal of Medicine and Science in Sports* 2008, 18:62-66.
- [20] Schween R, Gehring D, Gollhofer A: Immediate effects of an elastic knee sleeve on frontal plane gait biomechanics in knee osteoarthritis. *PLoS One* 2015, 10:e0115782.
- [21] Birmingham TB, Bryant DM, Giffin JR, Litchfield RB, Kramer JF, Donner A, Fowler PJ: A randomized controlled trial comparing the effectiveness of functional knee brace and neoprene sleeve use after anterior cruciate ligament reconstruction. *American Journal of Sports Medicine* 2008, 36:648-655.
- [22] Mohd Sharif NA, Goh S-L, Usman J, Wan Safwani WKZ: Biomechanical and functional efficacy of knee sleeves: A literature review. *Physical Therapy in Sport* 2017, 28:44-52.
- [23] Edwards PK, Ebert JR, Joss B, Ackland T, Annear P, Buelow J-U, Hewitt B: Patient Characteristics and Predictors of Return to Sport at 12 Months After Anterior Cruciate Ligament Reconstruction: The Importance of Patient Age and Postoperative Rehabilitation. *Orthopaedic Journal of Sports Medicine* 2018, 6:2325967118797575
- [24] Ebert JR, Webster KE, Edwards PK, Joss BK, D'Alessandro P, Janes G, Annear P: Current perspectives of Australian therapists on rehabilitation and return to sport after anterior cruciate ligament reconstruction: A survey. *Physical Therapy in Sport* 2019, 35:139-145.
- [25] Swearingen J, Lawrence E, Stevens J, Jackson C, Waggy C, Davis DS: Correlation of single leg vertical jump, single leg hop for distance, and single leg hop for time. *Physical Therapy in Sport* 2011, 12:194-198.

ACKNOWLEDGEMENT

We would like to thank the “School of Physiotherapy, University of Otago” for conducting the study.